
pekatherm
warm softness

- ES** INSTRUCCIONES DE USO
VENDAJES TÉRMICOS 2G
- EN** OPERATING INSTRUCTIONS
HEATING BANDAGES 2G
- RU** ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
БАНДАЖИ С ПОДОГРЕВОМ 2G
- PT** INSTRUÇÕES DE UTILIZAÇÃO
BANDAGENS AQUECEDORAS 2G

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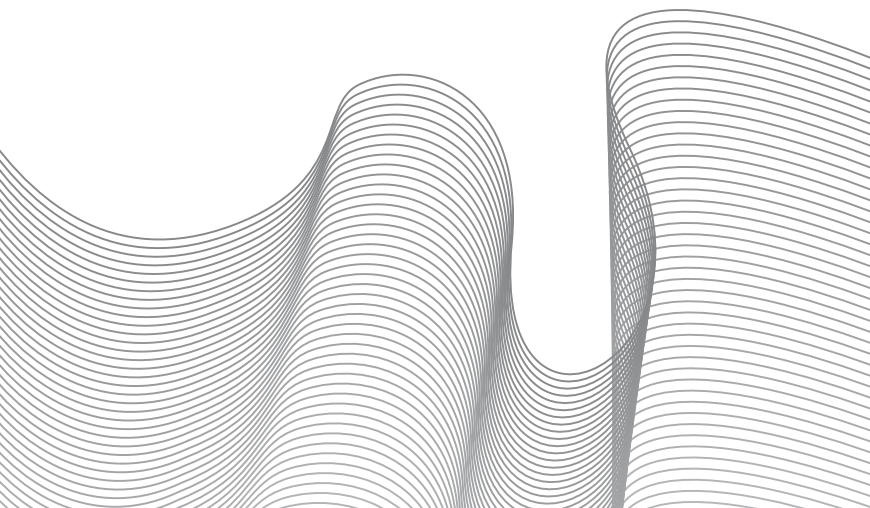IMPORTANT INFORMATION! PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE!

THESE INSTRUCTIONS PROVIDE INFORMATION AND ADVICE ON INSTALLATION, USE AND MAINTENANCE OF THE PRODUCT. READ THEM CAREFULLY BEFORE USING THE PRODUCT TO ENSURE SAFE AND PROPER USE. NON OBSERVANCE OF THESE INSTRUCTIONS CAN CAUSE PERSONAL INJURIES, AS WELL AS PRODUCT DAMAGE OR DESTRUCTION. KEEP THESE INSTRUCTIONS AND PASS THEM ON WITH THE PRODUCT IF SOLD OR LENT TO OTHERS.

PRODUCT PARTS

(see pic. 1, 2, 3, 4)

This set comprises heating bandage, battery, charger and user manual.

- a. Heating bandage (according to model)
- b. Battery
- c. Charger (plug may vary depending on the country)
- d. User manual
- e. Bandage connector

BATTERY PARTS

- f. Unlock button
- g. Red LED indicator
- h. Green LED indicator
- i. On/Off button
- j. Rear connector
- k. Charger assembly

BEFORE USE

1. After removing the packaging, inspect the product carefully to make sure it has no signs of damage caused during transport. In case of any doubt, do not use the product and take it to your dealer or to an authorized service centre.
2. Use only the battery and the charger supplied with the product. Never interchange them with those of other appliances.
3. This product must not be used with sick or disabled persons, with persons insensitive to heat or persons with restricted physical, sensory or mental capabilities. It also must not be used by children or by people lacking the required experience and knowledge of product operation unless they are properly instructed on how to operate the product safely. In all circumstances, while using the product, children should always be supervised by an adult to prevent them from playing with it.
4. Do not use the product under the effects of alcohol, sleeping pills or sedatives.
5. Do not use the product on animals.
6. Do not apply the product to the body areas where the skin presents scars, wounds, inflammation or any other alteration due to illness. The product should not be used by persons suffering from phlebitis, arthritis or other inflammations since heat may aggravate the disease. Consult your physician.
7. This product is designed for home use only. It is not intended for professional applications or medical use in hospitals.
8. Heating bandages should only be used for warming a specific part of human body. Any other use of the product will be considered hazardous. The manufacturer will decline all responsibility for any damage due to improper or incorrect use.

9. Pacemaker users should consult their physician before using the product.
10. Do not use the product folded or creased. Do not store the product creased.
11. Do not stick pins, needles or other sharp objects into the product.
12. Do not place heavy objects onto the product.
13. Never bring the product into contact with other sources of heat, such as other heating bandages, heating pads, heating underblankets, hot water bottles, etc.
14. Do not cover the battery while operating the product; do not place it on the heating element or under it while the unit is in operation.
15. Excessively long use of heating bandages running at a high temperature may produce skin burns. Avoid falling asleep with heating bandage in operation.
16. It is not recommended to apply the bandage directly on the skin for persons with sensitive skin. Take care to inspect frequently the skin in order to avoid skin burns. Disconnect the bandage immediately if you note burning or if the skin turns red or presents inflammations.
17. Do not leave the product unattended while in operation.
18. In any circumstances avoid wetting the battery and the charger. Do not use the product if wet.
19. Do not expose the product to environmental agents (rain, sun, etc.).
20. Keep the bandage, the battery and the charger away from fire.
21. Do not use creams or body oils that could damage the bandage.
22. Detach the battery and unplug the charger while not in use.
23. Detach the battery before cleaning the bandage.
24. For bandage cleaning please follow instructions in a specific section of this manual.
25. To detach the battery from the bandage or from the charger press the unlock button.
26. Do not hold the product by the battery to move or carry it.
27. Inspect the product frequently for signs of wear or damage. If these signs appear, do not use the product and take it to an authorized service centre.
28. This product has no parts intended to be repaired or replaced by the user. If any part of the product fails, or the product does not perform properly, or if there is any other kind of anomaly, do not attempt to repair it and instead take the product to an authorized service centre.
29. For any further information on our products please consult our customer service.

ADVICE ON BATTERY USE

30. Do not use the battery with other appliances.
31. Do not use the battery if the connector is damaged or presents signs of rust.
32. Never short-circuit the battery. Avoid keeping it together with small metallic objects. Do not wrap it in any conductor material.
33. Do not attempt to disassemble or alter any part of the battery.
34. Keep the battery away from children.
35. Keep the battery away from fire and water.
36. Battery is supplied half charged. Proceed to fully charge it before use.
37. To recharge the battery use only the supplied charger.
38. Once the battery is fully charged unplug it from the charger. Do not leave it plugged in when unnecessary.
39. Batteries gradually self-discharge with time. Proceed to fully charge the battery before each use.
40. Before storing the battery for a long period of time half charge it as explained in a specific section of this manual. Repeat the operation at least every half a year.
41. The capacity of rechargeable batteries gradually decreases. The supplied battery is designed to provide 300 discharge /charge cycles.
42. After a battery has reached the end of its lifetime take it to a green collecting point in order to assure the subsequent dispatch of the discarded elements for recycling, treatment and environmentally compatible disposal. Do not dispose it in fire or in water.

HEATING BANDAGES

Our heating bandages are developed and created in collaboration with physiotherapists and constitute the ultimate alternative in thermal therapy. It is a new line of high quality easy-to-use products which permit local targeted heat treatment. Low voltage operation assures maximum safety of use and its battery powered technology provides unlimited motion to the user. These products will actually make you forget the well known traditional heating products like hot water bottles, gel packs needing microwave or boiling water activating, etc.

- This new line of heating bandages includes seven models for different body parts:

MODEL	APPLICATION	MEASURES (CM)
AE802	Knee	25x19
AE804	Ankle	16x13
AE806	Lumbar	100x19
AE808	Shoulder	34x25
AE810	Elbow	25x14
AE812	Hand	23x13
AE814	Cervical	65x30

- All our bandages are uniquely designed to follow the complex shape of a specific body region providing an effective targeted heat treatment. All bandages are provided with built-in elastic straps with Velcro fasteners specially designed to secure the bandage onto a specific body joint in order to permit unlimited motion. Apart from the lumbar and the cervical models, all the bandages have universal anatomic design fitting equally for right and left member.
- An exclusive and specially designed rechargeable Li-ion battery is used to power all bandages. Ultra-lightweight construction of the battery reduces bulk and weight, while special connector provides an easy and quick lock. A modern built-in microprocessor controls the battery discharge/charge cycle, as well as the temperature of the heating element and overall running of the product, informing the user on bandage's operating conditions. Fast heating and temperature control performance are clearly superior to those of similar products on the market, yet the battery is utterly easy to use, with only one button, two LED indicators and audible signals. All models include a charger with replaceable plug.

- Three-layer cloth encloses the heating element and has a number of unique characteristics:
 - lightweight and fine, for increased wear comfort
 - softly textured thanks to its neoprene inner lining
 - provides high "next to skin" comfort thanks to its outer textile shell
 - highly elastic material permits perfect fitting and light compression on the joint where the bandage is applied
 - breathable material helps to maintain the skin dry during the treatment
 - insulating material retains heat and maintains the joint and adjacent muscles hot during treatment. Thermal insulation helps to prevent heat dissipation concentrating the heat in the affected area, which contributes to an easier heat penetration into the joint. Even if the bandage is used without battery the insulating material of the bandage retains body heat raising somewhat the temperature in the treated area in a natural way.
- Technical characteristics

Material	Neoprene
Power supply	Battery
Power	8W
Heating element	Resistive wire
Max. temperature, aprox.	55°C ^a
Battery autonomy	90min ^a
Fasteners	Elastic built-in straps with Velcro fasteners. Lumbar bandage is provided with a second pair of support straps
Battery characteristics	Li-ion 7.4V 650mAh. ON/OFF button, LED indicators, audible signals, fast lock
Charger	Wall adaptor Input voltage 90-240VAC 50-60Hz Output charge current 8.9V 325mA Fast charging Easy lock/unlock system
Charging time, aprox.	3h

^a Note: all measures were collected in specific laboratory conditions.

INSTRUCTIONS

BENEFITS OF HEAT THERAPY

Our products give heat, comfort and well-being easing muscle tension and stress. They also prove helpful in case of muscle or joint pains caused by an injury. Benefits of heat therapy are well known since ancient times. Local heat application stimulates metabolism and cell functions, dilates blood vessels increasing blood flow, oxygenation and nutrition of the affected area. This contributes to stimulate processes of regeneration and recovery of affected tissues, as well as elimination of toxins, products of tissue decay and lactic acid which improves vitality and helps to prevent other diseases. Heat reduces muscular tension and stiffness caused by overuse or overstretching and is also very helpful in case of rheumatic pains. Heat also relieves pain sensation since it reduces pain receptors stimulation and transmitting pain signals to the spinal cord and brain. Heat definitively makes great improvement in our sensation of well-being. Most athletes, professional and amateurs, perform some type of regular warm up before training and racing. A proper warm up can increase the blood flow to the working muscle which results in decreased muscle stiffness, less risk of injury and improved performance. In case of muscle injuries heat treatment is usually recommended after the initial inflammation is over.

REASONABLE HEAT TREATMENT

In order to bring positive effects, heat treatment does not require high temperatures or long time applications. It is with this premise in mind that the new series of heating bandages was created. Operating temperature is set at approximately 55°C, as measured in laboratory conditions. To get maximum benefit from heat treatment it is recommended to apply the bandage during 20 to 30 minutes. Consult your physician or physiotherapist.

There is wrong belief that the more heat and the more time the bandage is applied, the more will be the benefit. Heat must be clearly perceptible but never excessive. It should always be pleasant and should never cause physical trouble. In addition, in some circumstances heat may have a sedative effect. An excess of heat can lower blood pressure causing a sensation of dizzy or lightheaded or even faint.


If you perceive the heat emitted by your bandage is insufficient, try the following: remove the clothes from the part of the body to which the bandage is to be applied and then apply the bandage directly on the skin covering its other side with a towel. This will increase the heat transmission, as well as avoid heat dissipation. Take care to inspect frequently the skin in order to avoid skin burns. Disconnect the bandage immediately if you note burning or itch or if the skin turns red or presents inflammations.

USE INSTRUCTIONS

APPLICATION OF BANDAGES

AE802 (Knee)


Take the bandage with both hands and introduce your leg in it till the bandage's hole coincides with your knee cap. Battery connector should remain in the upper part of the bandage, above your knee.

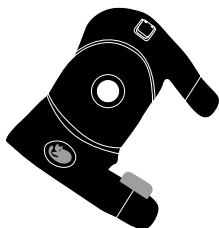Fasten the two elastic straps around your leg to fix the bandage.

AE804 (Ankle)

Take the bandage with both hands and put it on your foot making coincide the bandage's hole with your heel. Battery connector should remain at the rear part of the feet.

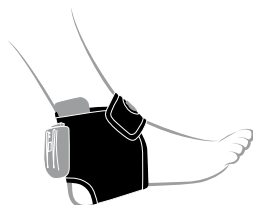
Fix the bandage with the help of the elastic strap.

INSTRUCTIONS

AE806 (Lumbar)

Put the bandage on your lower back making its stiffening rib coincides with your spinal column. Take the tags of the bandage and secure them onto your abdomen. Battery connector should be situated in the upper part of the bandage.

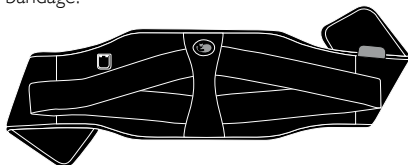Then take the elastic straps and pull them to your abdomen fixing them onto the bandage's surface. Thanks to its double fasteners the bandage will fit tightly on your lower back.

AE808 (Shoulder)

This bandage consists of two separate pieces: the proper bandage and the fixing strap.

1. Heating bandage

2. Fixing strap

Fix the strap to the bandage as shown below, depending on which shoulder the bandage is intended to be applied.

Introduce your arm in the bandage and place the bandage on your shoulder. Battery connector should be situated in the upper part of the bandage and the strap should remain on your back.

Right shoulder


Left shoulder

Take the loose end of the strap and pass it around your trunk fixing it on the front side of the bandage. Check that the bandage is fixed correctly. If necessary you can repeat the procedure adjusting the strap more tightly.

AE810 (Elbow)


Introduce your arm in the bandage until your elbow cap coincides with the bandage's hole. Battery connector should be situated in the upper part of the bandage, approximately above the triceps muscle.

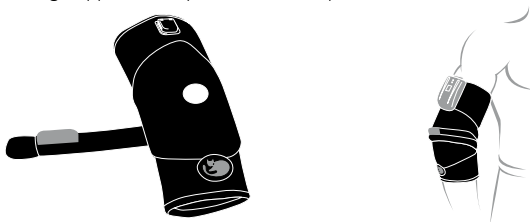Tighten the elastic strap around your arm and fix it with the help of Velcro fastener. Check that the bandage fits tightly on your elbow to ensure an effective treatment.

AE812 (Hand)


Introduce your hand in the the bandage, taking care to pass your thumb through the hole. Firmly fasten two straps with Velcro around your wrist.

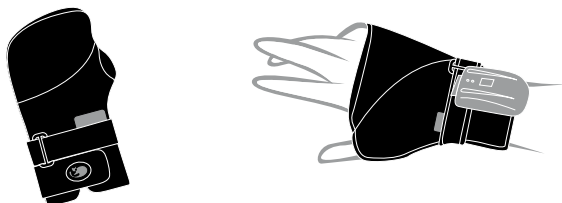AE814 (Cervical)

Take the bandage with both hands and place it on the back part of your neck, with the heating part of baddage in contact with your skin. Then take the strips at the ends of the bandage on either side of your neck and fasten the elastic strap with Velcro in the area of your throat.

Check that the heating part of bandage is properly centered on the back part of your neck, check that the contact between the bandage and your skin is good (without gaps) and check that the battery connector drops vertically over your chest. If needed, open again the elastic strap with Velcro and move the bandage to improve placement.

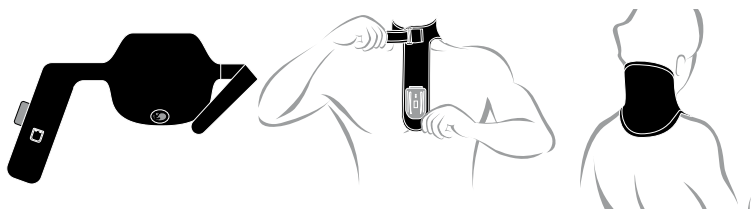INSTRUCTIONS

CONNECTING THE BATTERY

1. Take the battery and face carefully the connector on the back, with the connector on the bandage.
2. Press progressively the battery against the bandage connector till you hear an audible click.
3. Check that the battery has been properly connected

IMPORTANT NOTE: If the connection between the battery and the heating bandage is not correct, the bandage will not heat and the battery may drop to the ground and can be damaged.

To get a good connection, you must push hard enough the battery against the connector of the bandage, till you hear the audible click.

Check that the connection was successful, stretching slightly from the battery to be sure it is not easily disconnected from the bandage.

You also may find easier to connect the battery to the bandage before fixing the bandage on your body.

Use only the battery supplied with the product.

USE

1. Put the bandage onto a specific body part and fix it with the help of elastic straps (see above)
2. Connect the battery to the bandage, and check that the battery has been properly connected (see above)
3. Switch on the battery.
4. Green LED indicator will start flashing and a number of audio beeps will be emitted by the battery indicating the safety check is started. After the check is over the beeps stop and the green LED indicator remains lit. At this point the bandage begins warming up and in a few seconds you will be able to perceive a pleasant heat. Make sure the bandage is fixed correctly. Check that there are no creases or folds.
5. You may extend the treatment period if necessary; remember, however, that generally it should not last more than 20 to 30 minutes (see "Reasonable heat treatment"). Switch off the battery in any moment if you want to stop the treatment. Disconnect the battery from the bandage when not in use.
6. When the battery is completely discharged the control circuitry will switch off the bandage. Green LED indicator goes out and instead red indicator lights up and the battery emits audible signals. Disconnect the battery from the bandage and proceed to charge it.
7. After use let the bandage cool down before storage.
8. If the battery is not to be used for a considerable period of time store it disconnected from the bandage. Half charge the battery before storing as explained in a specific section of this manual.

BATTERY CHARGING

NOTE: Only use the battery supplied with the product.

1. Make sure your mains rating corresponds to the stated on the charger's label.
2. Switch off the battery and connect it to the charger pressing the battery against the charger connector till an audible click. Make sure the battery is connected correctly.
3. Plug the charger into the mains socket and switch on the battery. Red LED indicator will start flashing indicating the charging cycle is on.

NOTE: Solid green light indicates that the battery is fully charged. In this case unplug it from the charger. If no LED indicator is on or green LED is flashing slowly together with audio beeps, please check whether the battery is connected correctly to the charger and the charger is plugged correctly into the mains socket. Check as well if your mains socket has voltage. Switch off the battery and switch it on again after a few seconds.

4. Charging cycle will take about 3 hours for a completely discharged battery. When charging is over red LED indicator will stop flashing and green LED will remain lit. When the battery is fully charged switch it off and disconnect it from the charger. Unplug the charger from the mains to reduce consumption.

NOTE: Half charge the battery before storing if it is not to be used for a considerable period of time. Proceed as follows: completely discharge the battery connecting it to the bandage and then put it to charge for about 1 hour. Repeat the procedure at least once in half a year.

TECHNOLOGY pekatherm

The battery and its circuitry was specially developed and designed to power the new heating bandages to assure the best performance. At the start, when the circuit is switched on, the battery carries out an initial check recognizing the heating bandage and then proceeds to scan and store technical data which enables it to detect possible anomalies and to optimize the bandage's performance and assure the best operating conditions.

The battery is provided with high technology electronic circuitry, able to continuously control and regulate the temperature of the heating element. All the bandages share special fast heating technology.

TROUBLESHOOTING TABLE

PROBLEM	SOLUTION
1. After switching the battery on the green LED does not light and no audible beeps are emitted by the battery.	The battery may be completely discharged. Switch the battery off and turn it on after a few seconds. Repeat the procedure several times. Charge the battery.
2. After switching the battery on the green LED starts flashing and the battery emits audible beeps but after a few seconds the beeps stop, the LED indicator is off and the bandage does not warm up.	The battery may be connected to the bandage incorrectly. For safety reasons the battery will not operate if the connection is wrong. Check the connection. Detach the battery from the bandage and connect it again. You also may find easier to connect the battery to the bandage before fixing the bandage on your body.
3. After switching the battery on the red LED indicator is on and the battery possibly emits audible beeps.	The battery is partially or completely discharged. Disconnect the battery from the bandage and proceed to fully charge it as explained in this manual.
4. The bandage does not warm up.	Make sure the bandage is in operation checking the battery is switched on and green LED indicator is on. If necessary, consult numbers 1, 2 and 3 above. Remember that the bandage may have switched off automatically after the battery is completely discharged. Switch the battery off and then turn it on again. It may also be a question of feeling. Remember that quick checks carried out with uncovered bandage are not reliable. Try to cover the product with a blanket or a towel and wait at least 30 minutes. If you still do not perceive heat, take the product to an authorized service centre.

INSTRUCTIONS

5.The bandage warms up insufficiently.	If you perceive that the heat emitted by the product is insufficient, try the following: remove the clothes from the part of the body to which the bandage is intended to be applied and then apply the bandage directly on the skin, covering its other side with a towel. This will increase heat transmission, as well as avoid the dissipation of some of the heat. As it is explained in "Reasonable heat treatment" chapter local heat treatment does not require high temperatures.
6.After a period of proper functioning and before the battery is discharged the bandage seems to reduce notably its operating temperature.	With time skin gets used to higher temperatures and the feeling of heat may vary at the beginning of the treatment and some time after. Another possible reason of this performance may be the battery since the capacity of all rechargeable batteries gradually decreases. In no way it is a malfunction. If the battery capacity is the reason you may solve the problem providing a new battery for your bandage.
7.After a period of proper functioning green LED indicator turns red and the battery possibly emits audible beeps.	The battery has discharged. The bandage stops warming, red LED indicator lights up and the battery emits audible beeps. Proceed to charge the battery (See chapter on battery charging).
8.The battery appears to malfunction, LED indicators flash continuously and/or audible beeps are present.	The battery circuit may have found an anomaly and blocks the bandage operation for safety reasons. Switch off the battery, wait a few seconds and turn the battery on. If the problem persists take the product to an authorized service centre.
9.When the battery is put to charge green LED indicator lights up.	The battery is already charged. Disconnect it from the charger and unplug the charger from the mains socket.
10.When the battery is put to charge green LED starts flashing and the battery emits audible beeps. After a few seconds both LED indicator and audio beeps are off and the charging cycle fails to set up.	First check the connection between the battery and the charger. Switch the battery off, disconnect it from the charger and, after a few seconds connect it again. If the problem persists check if the charger is plugged in correctly. Check if your mains socket has voltage. If the problem persists take the product to an authorized service centre.

CLEANING

Before cleaning the bandage disconnect the battery. Never wet the battery and the charger.

For a basic cleaning of the bandages soft brushing is recommended. To get smaller stains out use a moist cloth or sponge with some liquid soap or cleaning foam. Apply the foam regularly on the product's surface and then brush it. If washing is necessary remember that only hand washing is permitted. Use a soft detergent and rub or brush softly avoiding wringing the bandage. Never bleach the bandage. After washing, extend the bandage flat or hang it, avoid folding it. Always place the bandage in the shadow since direct sunlight could damage the bandage. Do not iron the bandage.

STORAGE

Let the product cool down before folding it for storage. Avoid excessive folds or creases. Store it in a cool dry place. Do not place any heavy objects on the product. Follow the instructions on battery storage when not in use for long periods of time. Inspect the product (the bandage, the battery and the charger) frequently for signs of wear, damage or misuse. If these signs appear, do not use the product. Remember this product has no parts intended to be repaired or replaced by the user. Do not try to carry out repairs and, instead, take the product to an authorized service centre.

Repairs may only be carried out by our after sales service or at an authorized service centre. Improper repair may cause considerable danger to the user. In the event of any doubt, please contact our customer service.

ADVICE ON DISPOSAL

According to EU Directive 2002/96/EC WEEE (Waste of Electrical and Electronic Equipment), electrical and electronic devices must be disposed of separately from household refuse, delivering them to a special refuse collection point, in order to assure the subsequent dispatch of the discarded equipment for recycling, treatment and environmentally compatible disposal, thus contributing to the prevention of possible negative effects on the environment and on health.

Your local authorities will provide you further information on how to dispose of your inoperative appliances.

WARRANTY

pekatherm S.L. guarantees this product against manufacturing and material defects for a period of two years from the date of purchase.

pekatherm S.L. will replace or repair the defective item free of charge, including labour and spare parts, if the defect derives from manufacturing and/or material reasons, and if the product is handed by the customer to an authorized service centre.

If replacement or repair is not possible, the customer may claim for refund or buy a superior model paying the difference.

This warranty will not be valid if not accompanied by the original purchase bill or ticket. Application of the warranty will not imply any extension of its period.

This warranty will not cover any defect due to reasons beyond manufacturing, such as transportation damages, non observance of the operating instructions, careless or improper use, wear and tear, damages caused by a blow or fall, water, chemical agents, etc.

This warranty will not cover any cases of improper use of the product, such as professional or medical use.

This warranty will not cover any damage if improper repairs were attempted by unauthorized personnel.

pekatherm S.L. declines all responsibility for any damage caused directly or indirectly to persons, pets or properties deriving from non observance of these operating instructions.

This warranty is only valid in the country where the product was purchased.

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