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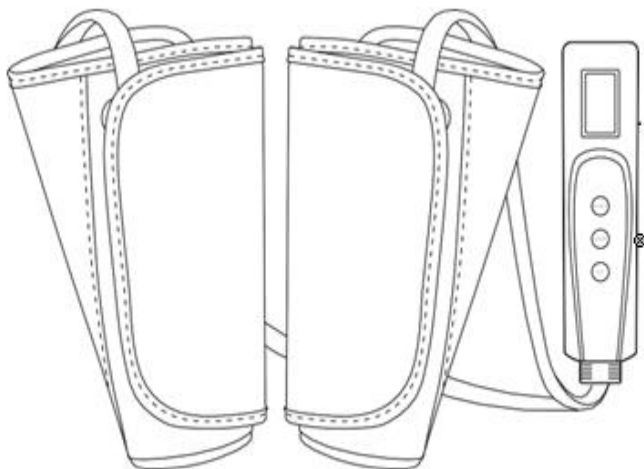
Area of use

Massager MFC-50 is used for the procedure of pressotherapy in the area of the thighs, calf muscles, ankles, feet, hands and forearms. The massager is specifically designed for easy and convenient procedures at home. A powerful compressor pumps alternately air in the air chambers along the length of the cuffs. Due to the Velcro fasteners on the cuffs tight wrapping and fit is provided. Compression massage improves blood circulation, lymph outflow, helps to relieve muscle tension and frees the limbs from edematous fluid. This massager promotes weight loss and it is effective in the fight against cellulite. The function of heating stimulates blood circulation. MFC-50 is ideal for people who are long at work on their feet or have discomfort in feet, for people who lead a sedentary lifestyle.

The massager is intended for personal use. It is not intended for commercial use of any kind, nor for use in medical institutions. To use all the advantages of your massager carefully study this user manual and follow the rules of care and safety. This will extend the service life of your massager.

Safety instructions

Before using of the massager read the user manual and safety instructions carefully.

- Ensure that all parts of the device are in good condition before use. If a defect or fault is revealed, turn off the device. Only specialists of authorized services can repair the device. Never open the device yourself. Separate components of this device cannot be disassembled, repaired or reinstalled.
- Use the device only in absolutely dry rooms, making sure that your body is also dry. It is forbidden to use the device in the bathroom, shower, above sink filled with water, in swimming pools, sauna and outdoors.
- Never use the device in close proximity to heat sources.
- Do not use the device more than 20 minutes.
- This device is forbidden to be used by persons (especially children) with physical, sensory or mental disabilities, as well as persons who do not have sufficient experience or knowledge to handle the device, except when using the device under the control of the responsible person or after obtaining the necessary information.
- Ensure that children do not use the device as a toy.
- Always unplug the device before cleaning. During cleaning make sure that liquids do not hit the switches and buttons. The device can only be used when it is completely dry.
- It is forbidden to disassemble the device. It has no parts that the user can repair.
- It is forbidden to remove the panels of the device, since the repair of internal components can only be carried out by qualified specialists.

- Do not insert needles or other sharp objects into the massager.
- Never use this device in conjunction with electro-medical equipment, electronic life support systems such as ventilators or portable medical devices such as an electrocardiograph.
This device can cause disturbances in the operation of the above-mentioned medical equipment.
- Do not leave the working device under pillows, blankets, etc.
- If the power cord is damaged, it must be replaced in the service center or by a qualified person to avoid danger.

Warning

Consult your doctor before using the device if :

- You are pregnant;
- You use a pacemaker and others;
- You suffer from one of the listed diseases: influenza, fever, varicose veins, thrombosis, hepatitis, diabetes, inflammation, ischias, loss of intervertebral discs, etc.;
- You undergo rehabilitation after surgery.

Do not use the device if there are open wounds, burns, hematomas, tumors, etc. on your body.

Do not use the device when under the influence of alcohol.

Stop massage if you have unpleasant or painful feelings.

If you feel limb pain during or after using the device, immediately stop using it and visit your doctor.

The device has heating surfaces. Persons who are sensitive to the effects of heat should take into account that the useful life of the device is 20 minutes and the device can only be used under supervision.

Care

Turn off the device and let it cool down before cleaning.

Use cloth wetted with soft detergent for cleaning.

The device can be used again after it has completely dried.

It is forbidden to remove the device panels or open fasteners because the repair of internal components can be carried out only by qualified specialists.


Attention!

Never use abrasive or aggressive cleaning agents and do not put the device into the water.

Description of the device

Box contents:

1. Massager
2. Adapter
3. User Manual with warranty card

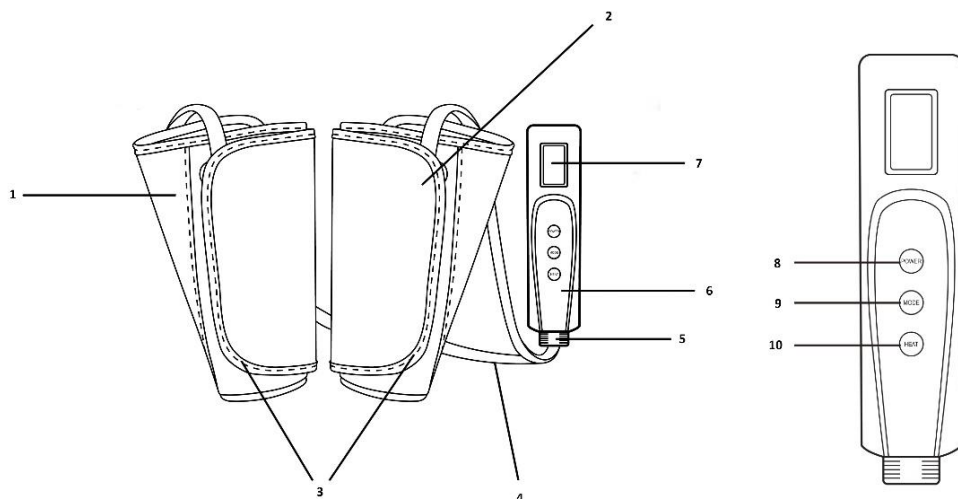1. Right cuff R
2. Left cuff L
3. Velcro fasteners
4. Air pressure hose
5. Place of connection of hoses
6. Control panel
7. LCD screen

Description of control panel :

8. ON/OFF and intensity choosing button(POWER)
9. Mode choosing button (MODE)
10. ON/OFF of heat button (HEAT)

Operating the appliance

1. Make sure that the device, plug and cable are not damaged
2. Put the cuffs on your arms or legs, covering the calf muscles, thighs or feet, as shown in the pictures below:

3. Get into a comfortable position. Fix the massager on your legs or hands with the help of Velcro fasteners.
4. Connect the massager according to the pictures: the hose from the left cuff L to the L socket on the control panel and the hose from the right R cuff to the R socket on the control panel.
5. Connect the device to the network.

6. Press the POWER button for 3 seconds. The massager will start working. By default, the first intensity level will be set - the lowest. Pressing again will set the intensity to medium, pressing third time - to high intensity.
7. Press the MODE button to choose the massage mode. MODE 1 will be set by default. Pressing the MODE button again, you can choose one of the three modes (the LCD display will show the mode you choose):

MODE 1-compression in the upper and lower part of the cuff

MODE 2-compression in the upper part of the cuff

MODE 3-compression in the lower part of the cuff

8. Additionally you can turn on heating by pressing the HEAT button . When you press the HEAT button, the first heating level will be set – the lowest. The second time you press it, the medium level will be set, and the third time you press it, the high level of heating will be set. When the HEAT button is pressed for the fourth time, the heating mode will be switched off.

9. The device will automatically switch off after 20 minutes of continuous working.

10.To switch off the device before automatic shutdown time, press and hold the POWER button for 3 seconds. The compressor will stop working, the pressure in the cuffs will be lowered. After that, the cuffs can be removed.

Wait until the device completely cools down before using it again.

It is not recommended to use the massage device for more than 20 minutes.

Attention!

- If a problem is revealed, stop using the device immediately and forward it to the nearest service center for testing, replacement or repair.
- For safety reasons all repair and maintenance of the device must be carried out only by representatives of authorized service centers.

Do not try to fix the problem yourself!

Possible problems and how to fix them

Problem	Possible cause	How to fix
Device does not switch on.	The device is not connected to a power source.	Make sure that the device is connected to power.
	The device was switched on not correctly.	Press the POWER button for 3 seconds
The device suddenly stopped working.	You use the massager more than 20 minutes.	It is normal. The device works properly.
During using you hear sounds .	It can be a sound of air compressor.	
The device works not correctly.	The device was kept at low temperature for some time.	Leave the massager for 2 hours at room temperature, then use it.
The massager works but pressure is low.	The chosen level of intensity is too low for you.	Press the button of intensity and change the level of intensity.
	There is an air leak.	Make sure that both hoses are tightly connected to the control panel.
The sound of escaping air is heard.	There is correction of the pressure in the air chambers.	It is normal situation for working of the massager.

Technical specification


Model: MFC-50

Rated voltage: AC 100V-240V / 50-60Hz

DC 12V/1.5 A

Power: 18W

Cuff Dimensions:

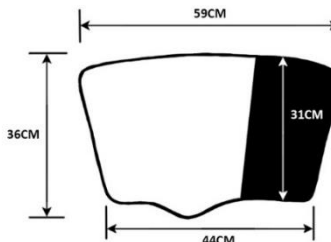Air pressure:

Pressure in the upper part of the cuff –

at low intensity - 28 kPa (0.21 mm Hg)

at medium intensity-44 kPa(0.33 mm Hg)

at high intensity-48 kPa(0.36 mm Hg)

Pressure in the lower part of the cuff –

at low intensity - 42 kPa (0.32 mm Hg)

at medium intensity-58 kPa (0.44 mm Hg)

at high intensity-62 kPa(0.47 mm Hg)

Weight: 0.8 kg

Device material:

Control panel-100% ABS plastic

Air pressure hose -100% silicone

Cuff- 100% polyester

The manufacturer reserves the right to change parts of the product without prior notice for different schedule lines without affecting the basic technical parameters of the product. This can entail a change in the weight and dimensions of the product, but not more than +/- 5-10%.

Disposal of electrical equipment

In the interests of the environment, after the end of the service life, the device cannot be disposed of together with household waste. Disposal shall be carried out through appropriate collection and processing points for electrical and electronic equipment. If you have questions, contact the local municipal service responsible for waste disposal.

Warranty card

Warranty Terms

- 1.If a fault is revealed during the Warranty Period, the Product shall be exchanged with the retailer that sold the Product.
2. Warranty replacement shall be carried out if the warranty card is filled in and subject to the operating rules described in the instruction manual.
- 3.Product warranty does not apply in the following cases:
 - mechanical damages;
 - failure of the article due to ingress of foreign objects and liquids, insects, etc.;
 - use of the article in conditions and modes other than domestic.
- 4.The warranty does not apply to accessories and components.
5. The warranty also loses force if during the warranty period the repair of the faulty product was carried out by non-authorized persons.
- 6.The warranty period is 1 year.
- 7.The service life is 3 years.

The date of production is indicated on the package.

Manufacturer: Xiamen Healthpal Electronic Co.,Ltd.

Address: No. 148 Binglang Xili, Siming District, Xiamen, Fujian, China

Importer/Authorized Organization: LLC "PALLADA" 111675,Russia, Moscow, Rudnevka Street, House 4, P IV K 1.2 OF 1 RM 2/5, Russia

Customer service support:
<https://www.planta.ru.com>

+7 (495) 137-80-32

Product is certified

Warranty card

Device: _____

Model: _____

Date of sale: _____

Selling company information: _____

I received the product, I have no complaints about the configuration and appearance, I agree with the terms of the guarantee. Buyer's full name and signature: : _____

Opened the package, checked and sold. Seller's name: _____

Stamp