# A E N O

User Manual


## **Toothbrush** DB5/DB6


(Model: ADB0005 / ADB0006)

29.11.2021 | v 1.0




EN.....	6
RU.....	11
AM.....	27
AR.....	22
BG.....	26
CN.....	33
CZ.....	36
DE.....	41
EE.....	47
FR.....	52
GE.....	58
HU.....	64
KZ.....	69
LT.....	75
LV.....	80
PL.....	85
RO.....	91
SK.....	96
UA.....	101
UZ.....	106

a)

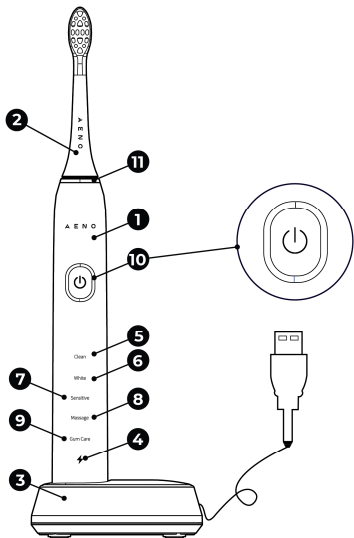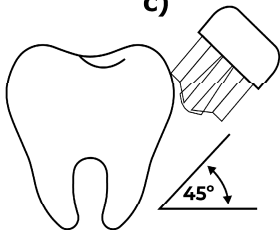**b)****c)**

d)

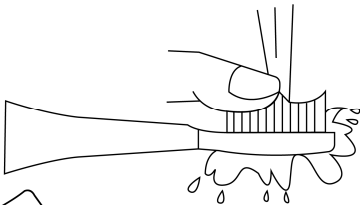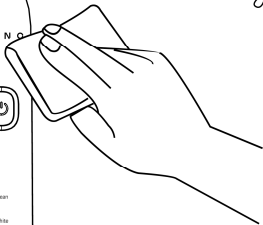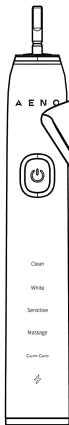**EN** Toothbrush **AENO DB5 / DB6** is designed for daily cleaning of teeth and gums.

### Specifications

**Toothbrush:** non-removable lithium-ion battery (DC) 3.7 V, 750 mAh. Power 2.0 W. Noise level (0.3m distance): ≤65 dB. Time to fully charge the battery 4 h ± 30 min. Working time after full charge: not less than 150 min. The degree of protection of the enclosure is IPX7. Bristle hardness: medium. LED indicators for operating modes: Clean, White, Sensitive, Massage, Gum Care. LED indicators battery status: charged, discharged. Automatic timing function. The degree of protection of the enclosure is IPX7. Dimensions of the brush with attachment (D×H): Ø28 mm × 246 mm.

**Charging station:** Input voltage: (DC) 5.0 V, 1.0 A. Power 3.0 W. The degree of protection of the enclosure is IPX7. Dimensions (L×W×H): 80mm × 48mm × 26mm.

### Scope of delivery (fig. A)

Toothbrush **AENO DB5** (white) / **DB6** (black) (handle (1), 2 replaceable toothbrush heads (2)), charging station (3) with a USB-A power cord (1.2 m), a quick user's guide, a warranty card.

### Limitations and Warnings

The appliance is not intended for use by children or persons with reduced physical/mental/intellectual capabilities unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Do not leave a plugged-in device unattended. Do not turn on the damaged device. Disconnect the device from the mains by holding the power plug, not the cord; do not use an extension cord. The toothbrush is designed exclusively for cleaning teeth and gums. After surgical intervention, such as surgery on the gums or in the oral cavity, as well as in the presence of gum and oral cavity diseases, long-term dental problems or suspicions of similar problems (fillings, crowns, etc.), you must first consult your

dentist about the possibility of using the device. If your gums begin to bleed after using the device and the bleeding continues for more than seven (7) days, contact your dentist immediately. Charge your toothbrush only on the original charging cradle, placed on a level and stable surface. If the charging cradle accidentally falls into water, you must first unplug it from the mains and only then remove it from the water and dry it thoroughly. The area where the mains cable of the charging cradle and the USB cable are located must be dry. Each family member should use their own brush head. Replace the brush head at least once every three months, depending on the degree of wear. Use only original replacement tips.

**ATTENTION!** Do not wash the device or any part of it in the dishwasher.

### **Preparing the device for use**

Unpack the device and place the charging station (3) on a stable horizontal surface. Connect the charging station to the mains using a USB cable and an AC adapter (not supplied). Place the device (1) on the charging station. The charging LED (4) will turn red. When the LED indicator (4) glows green, it means that the battery is fully charged, and the device is ready for use.

**ATTENTION!** The AC adapter is not included in the package. You can use the AC adapter from your phone, tablet, or other device, if it provides power options 5.0 V (DC), 1.0 A.

To prolong battery life, charge the brush for at least 20 hours before first use. Wait until the battery is completely discharged before the second charge. The second charge should last at least 6 hours.

At the end of charging, you must unplug the charging station from the outlet. When the battery is fully charged, the brush runs for about 150 minutes. If you brush your teeth twice a day for 2 minutes, the device will work without charging for up to 1 month. When the battery is discharged, the charge indicator (4) flashes red and starts flashing slowly.

**NOTE.** When charging for the first time or after a long interruption (about 6 months), the LED indicator (4) may turn on in about 10 minutes after starting charging.

### Operation of the device

Place the cleaning pad (2) onto the handle of the device (1) as shown in figure **b**. Apply some dentifrice to the brush bristles. Press the power button (10). One of the indicators starts to light, showing the mode of operation of the device. The device starts working in this mode. To turn off the device, press the button (10) again.

Every 30 sec. a vibration signal is triggered, reminding you of the opportunity to change the operating mode. To change the mode, switch the device off and on again with the button (10). Then quickly press the button (10) several times in succession, selecting the desired mode. Each time you press, the mode changes to the next one. After each use, the toothbrush remembers the last mode. This mode will be automatically selected the next time the device is turned on.

Guide the brush head at a 45° angle to the surface of the teeth, as shown in Figure **c**. In this position, the toothbrush head works most efficiently.

### Device operating modes table

Индикатор (рисунок <b>a</b> )	Режим	Описание
<b>5</b>	«Clean»	Daily standard cleaning mode
<b>6</b>	«White»	Teeth polishing and whitening mode. Recommended for smokers and frequent coffee drinkers
<b>7</b>	«Sensitive»	Mode for delicate, gentle cleaning of gums and teeth

8	«Massage»	Intensive gum massage. Recommended for healthy gums after standard cleaning
9	«Gum Care»	Gentle massage of the gums. Recommended for sensitive gums after standard cleaning

### Auto power off mode

After 2 minutes of operation, the device automatically turns off. You can turn on the device again by pressing the button (10).

### Maintenance

Replaceable toothbrush heads should be washed in running water at a comfortable temperature for you, as shown in Figure **d**. Try not to completely submerge the toothbrush head in water. Do not use aggressive or cleaning agents when washing the attachment. It is allowed to treat the bristles of the toothbrush head with a mouthwash from time to time. Any dirt on the handle of the device and the charging station should be removed with a damp cloth, as shown in figure **d**. Then you should carefully wipe the surface with a dry cloth. For stubborn dirt, apply dishwashing liquid to a cloth to clean the handle of the device. Remove any remaining moisture with a towel after washing.

**NOTE.** Perform maintenance on the device and its components only after turning off the device and disconnecting the charging base from the mains.

### Elimination of possible malfunctions

- The toothbrush does not turn on. Possible cause: The device battery is low. Solution: charge the device for at least 4 hours.
- The device placed in the charging station does not charge. Possible reasons: damage to the cord of the charging base, incorrect parameters of the power adapter. Solution: Make sure the cord is intact and that the outlet is working. Check if the AC adapter has (DC) 5.0 V, 1.0 A output.

- The brush head vibrates very weakly. Possible Cause: The attachment is too tight on the handle of the device. Solution to the problem: make sure that there is a small gap at the mating point (11), about 1 mm.
- Short working time after fully charged. Possible cause: The device has not been used for more than 6 months, so the battery may not hold a charge well. Solution: contact the service center and replace the device battery.

**ATTENTION!** If none of the possible remedies resolves your problem, contact your supplier or service center. Do not disassemble the device or try to repair it yourself.

### Recycling information

These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment.

To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility.

You should contact your local household waste disposal service for details.

The warranty period and service life shall be 2 years as from the date of purchase of the device.

Manufacturer: ASBISc Enterprises PLC, 43 Kolonakiou Street, Diamond Court, Ayios Athanasios, 4103 Limassol, Cyprus. Made in China. For current information and details on the device description and specification, as well as connection process, certificates, warranty and quality issues, see relevant Installation and Operation Manuals available for downloading at [aeno.com/documents](http://aeno.com/documents). All trademarks and names herein are the property of their respective owners

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