
PLANTA

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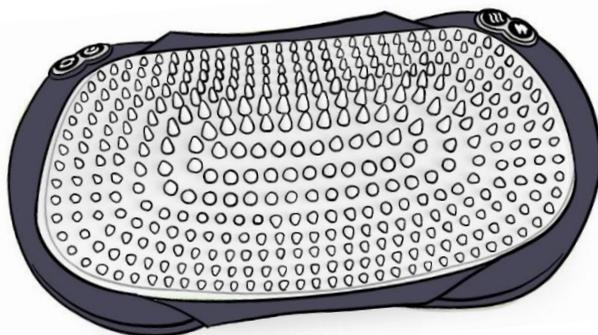MANUAL

ATTENTION! SAVE THIS MANUAL DURING ALL USING!

With PLANTA MP-021 Acupuncture massage pillow you can enjoy massage at anytime and anywhere. Massager can be used to massage any area of the body - back, shoulders, lower back, legs, hips, abdomen. This massage pad has a wide range of applications. Massage is one of the most effective practices for improving the body. Mankind does not know the more ancient way of getting rid of ailments - thousands of years ago, people treated each other. Usually, massage serves to relax muscles and improve well-being. Massage with this device has a relaxing and stimulating effect. It has been proven to reduce muscle strain and relieve fatigue. Thanks to acupuncture, massage becomes more effective and high-quality, even better working out trigger points, muscles of the whole body. Device is only for personal use. It is not intended for commercial use of any kind, nor for use in medical institutions. To take full advantage of your massage pillow, carefully review the operating manual and follow the rules of care and safety. This will extend the service life of your massager. Save this operating guide for future use or when transferring the device to another person. For more information, visit: www.planta.ru.com

SAFETY REQUIRMENTS

* Before using the device, consult your doctor if you:

- pregnant;
- use a pacemaker, etc.;
- suffer from one of the following diseases: flu, fever, varicose veins, thrombosis, hepatitis diabetes, inflammation, sciatica, intervertebral disc prolapse, etc.;
- undergo rehabilitation after surgery.

* Do not use the device if there are open wounds, burns, bruises on the body, tumors, etc.

* In addition, you can not use the device to treat diseases of the joints, spine, head and other sensitive areas.

* Do not use the device while under the influence of alcohol.

* Stop the massage if you feel unpleasant or painful sensations.

* If you have muscle or joint pain after using the device, stop using it immediately and consult a doctor.

* Do not leave the device running under pillows, blankets, etc.

* If the power cord is damaged, it must be replaced at a service center or by a qualified technician to avoid danger.

* The device has heating surfaces. Persons sensitive to heat should take into account that the useful life of the device is 45 minutes and the device can only be used under supervision.

Notes

*If you use the device just before going to bed, remember that massage it can have a stimulating effect and cause insomnia.


*The massager turns off automatically after 15 minutes, which corresponds to the optimal the duration of the massage, which has a beneficial effect on the body.

Care


- Use a cloth moistened with a mild detergent for cleaning.
- Possible risks
- Before cleaning, pull the plug out of the socket and let the device cool down.
- The device can be used again after it is completely)
- it will dry out.
- It is forbidden to remove the panels of the device or open the fasteners,
- because the repair of internal components can only be carried
- out by qualified specialists.

Attention! In no case do not use abrasive or aggressive cleaning agents and do not immerse the device in water.

DEVICE DETAILS

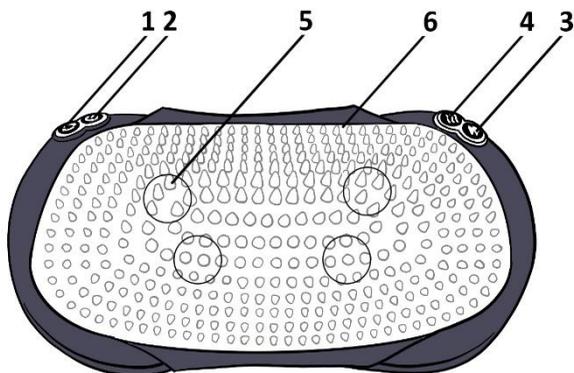2 POWER ADAPTERS

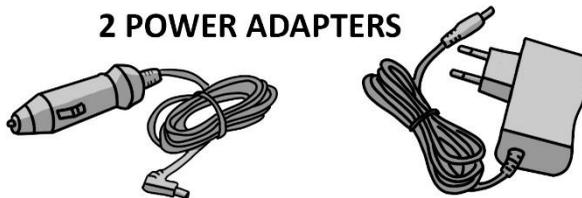
1. Enable button (On/Off)
2. Button for selection of direction of massage rollers rotation (clockwise or counterclockwise)
3. Button for selection of rotation speed of massage rollers (3 intensity options)
4. Heating. Button enables or disables heating in the area of massage rollers
5. Massage rollers (Shiatsu)
6. Acupuncture cap. The cap is attached to the massage cushion using a lipstick (see Fig below)

The massager can be used with this cloak to enhance the effect of massage and acupuncture, or without it to obtain a classic Shiatsu massage.

POWER ADAPTERS

This device can be used both in the car and in the home. The kit has 2 adapters: for the car (12V), for working from the network (220V).

To connect, connect the device plug to the selected power adapter. (see. Rice. below)

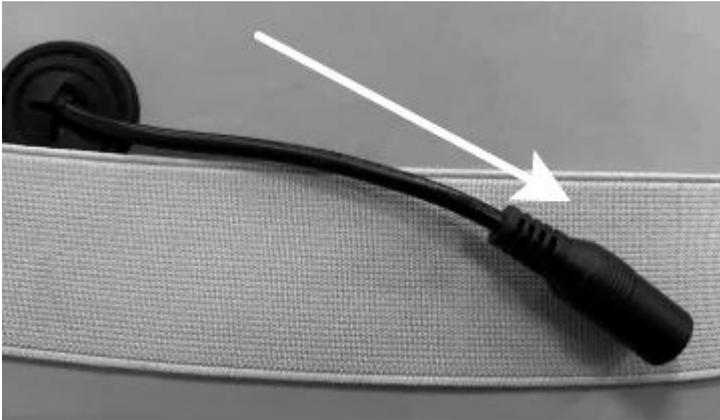USING


Make sure that the instrument, plug and cable are intact.

Place the massage device on the chair or place it conveniently on the car seat.

Connect the device to the network.

Attention! Acupuncture cap on sticking stickers, and easily removed, allowing you to use the device separately from the accessory at your discretion.

1. First, select the type of massage (with acupuncture capping/without acupuncture capping), then carefully press the area you want to massage to the massage device.


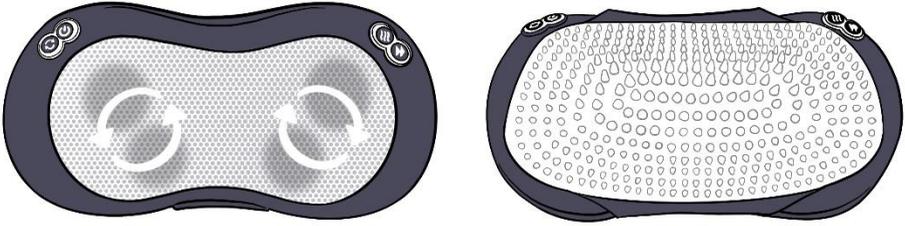
Attention! The acupuncture nozzle can be used separately from the massager.

2. The device can be used on the hips, back, hands, feet, lower back and neck-collar zone. Make sure that the massager position is convenient for you.

Massage should always be relaxing and pleasant. If the massage causes unpleasant sensations or pain, interrupt the massage, change the position or force of pressure.

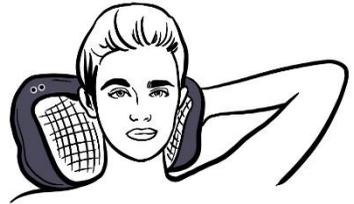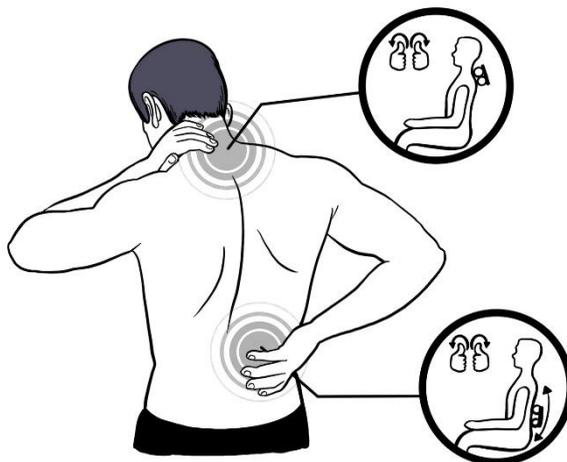3. To turn on the massage device, press the power button. Massage heads will begin to rotate. To change the massage intensity, click the speed switch button. By default, the average speed will be set. Pressing the button again will slow down the rotation. With the third press, the massage speed will increase.
4. Press the heat button to activate the heat function. Pressing the button again will disable this function.
5. Change the direction of rotation of the massage rollers using the rotation direction switch.
6. Use a massage device for no longer than 15 minutes. Excessive muscle stimulation during prolonged massage can lead to their tension instead of relaxation.
6. When the massage is finished, the device will automatically turn off (in 15 minutes). To abort the session, press the power button again and the device will disconnect.

STORAGE, TRANSPORTATION AND REGULAR INSPECTION

- If you do not use the massager, store it in a dry place.
- Before folding the storage device, let it cool completely.
- Do not place heavy objects on the device during storage, otherwise parts of the device may break.
- Do not allow the device to fall. This may cause damage to it.
- Be careful when transporting the massage cape.

Faults

Problem	Reason	Solution
The device does not turn on	The device is not connected to the power supply	Check the correct connection of the massager
	There is no voltage in the network	Connect the massager to a known serviceable outlet
Massage rollers move too slowly	The level of massage intensity is too low	Increase the intensity of the massage
	Too much pressure on the massage nodes	Reduce the pressure on the massage nodes
The massager is not working properly	he device has been at a low temperature for some time	Leave the massager at room temperature for 2 hours, then turn it on again

TECHNICAL SPECIFICATIONS

Model: MP-021

Power: 100V-240V/50-60Hz

DC 12V/2A


Power: 24W

Pillow size: 18 * 12 * 35 cm

Weight: 1.28 kg

Cap acupuncture size: 15 * 10 cm

Number of spikes: 350 pcs

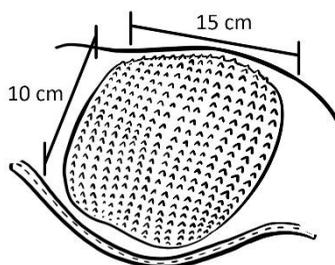The manufacturer reserves the right for different delivery lots without prior notice to change the component parts of the product, without affecting the main technical parameters of the product. This may entail a change in the weight and dimensions of the product, but not more than + \ - 5-10%.

Disposal

In the interest of environmental protection, at the end of its useful life, the appliance must not be disposed of with household waste. Disposal must be carried out at the applicable collection point for the recycling of electrical and electronic equipment. If you have any questions, please contact your local waste disposal company.

WARRANTY CARD

WARRANTY TERMS

1. If a fault is detected during the Warranty Period, the Product shall be exchanged with the retailer that sold the Product.

2. Warranty replacement shall be carried out if the warranty ticket is filled in and subject to the operating rules described in the instruction manual.

3. Product warranty does not apply in the following cases: - mechanical damages; - failure of the article due to ingress of foreign objects and liquids, insects, etc.; - use of the article in conditions and modes other than domestic.

4. The warranty does not apply to accessories and components.

5. The warranty also loses force if during the warranty period the repair of the faulty product was carried out by non-authorized persons.

6. The warranty period is 1 year.

7. The service life is 3 years.

The date of manufacture is indicated on the package..
 Manufacturer: Xiamen Healthpal Electronic Co.,Ltd.

Address: No. 148 Binglang Xili, Siming District, Xiamen, Fujian, China

Importer/Authorized Organization:
 LLC "PALLADA" 111675 MOSCOW,
 RUDNEVKA STREET, HOUSE 4,
 P IV K 1.2 OF 1 RM 2/5, RUSSIA

Customer Service: <https://www.planta.ru.com>

+7 (495) 137-80-32

The goods certified

WARRANTY CARD

Product:

Model:

Date of sale:

Selling company information:

Product received, claims against configuration and appearance I do not mean, I agree with the terms of the guarantee. F.I.O. and buyer's signature:

Opened the package, checked it and sold it. F.I.O. of the seller:

Stamp