
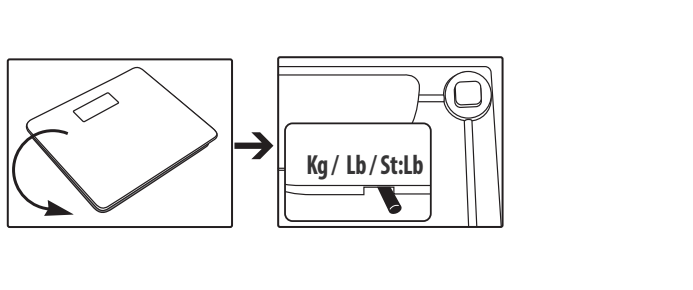
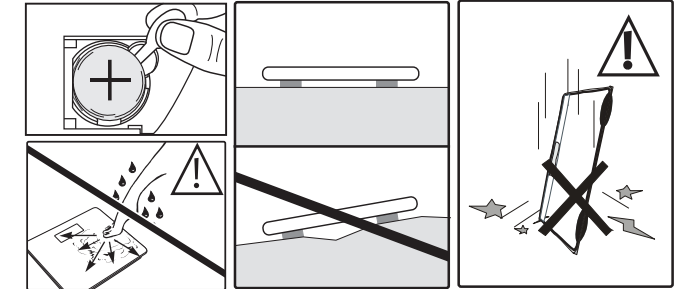
Aprilla

EN Instruction Manual
FR Mode D'emploi
NL Haarborstel
DE Bedienungsanleitung
IT Manuale Utente
ES Manual De Usuario
PT Manual De Utilizador

TR Kullanım Kilavuzu
HR Korisnički Priručnik
RS Uputstvo Za Korištenje
SI Navodila Za Uporabo
RU Инструкция
KZ нұсқаулық
EN Bathroom Scale
FR Pèse-personne
NL Weegschaal
DE Personenwaage
IT Bilancia pesapersona
ES Báscula de baño
PT Balança de banheiro
TR Banyo Baskülü
HR Osobna Vaga
RS Kupaonske Vage
SI Vaga Za Kupatilo
RU Весы Для Ванной
KZ Ванна Таразы

Parts Description / Description Des Pièces / Teilebeschreibung / Descrizione Delle Parti / Descripción De Las Piezas / Descrição Dos Componentes / Parça Tanımları / Описание частей / Opis Dijelova / Opis Delov / أجزاء الوصف

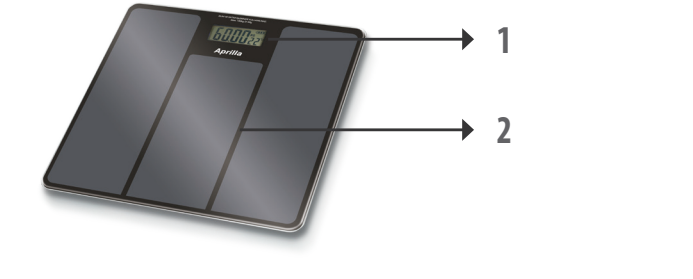Technical Data / Données Techniques Technische Daten / Dati Tecnici / Datos Técnicos / Dados Técnicos / Teknik Bilgi / Технические данные / Техникалық деректер

Power supply / Source de courant / Источник питания / Нәр беруші / Energieversorgung / Alimentazione elettrica / Fuente de alimentación / Fonte de energia / Güç kaynağı / Napajanje

CR2032/3V

Capacity / Capacité / Capaciteit / Kapazität / Capacità / Capacidad / Capacidade / Kapasite / Kapacitet / Емкость / Сыйымдылығы
180 kg max. / D: 50 g

Aprilla-Care.com

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EN Instruction Manual
SAFETY
• The manufacturer cannot be held responsible for damage if the safety instructions are not observed.
• This appliance is only to be used for household purposes and only for the purpose it is made for.

DEVICE SAFETY:
• The appliance must be placed on a stable, level surface.
• Do not step with wet feet. Its surface is slippery. It may cause injury.
• Do not drop the device.

BATTERY WARNINGS:
• Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
• Keep batteries away from children. They are not a toy!
• Do not open the batteries by force.
• Avoid contact with metallic objects.(Rings, nails, screws et cetera)there is a danger of short-circuiting!
• As a result of a short-circuit batteries may heat up consider-ably or even catch fire.this may result in burns.
• For your safety the battery poles should be covered with adhesive strips during transport.
• If a battery leaks, do not rub your eyes with your hands. If in contact with the liquid, the hands should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.
• Attention should be drawn to environmental aspects of battery disposal.

PARTS DESCRIPTION
1. Digital screen 2. Glass Surface

BEFORE FIRST USE
• Take the appliance and accessories out the box. Remove the stickers, protective foil or plastic from the device.
• Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
• It is advisable to take measurement at the same time of the day and at the same place.
• Step on the weighing platform gently and when you are measuring, please keep still.

USE
• The maximum weight capacity is 180kg. The minimum measurable weight is 3kg. Graduation is 50gr.
• Select with the little knob aside the battery holder the weight measuring unit KG, LB or ST:LB.
• Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate measurements.
• Stand on the scales with both feet. The device switches on automatically.
• Within seconds, personal weight is established electronically and displayed.
• Now for example you can take up your baby, if the weight is heavier than 3 kg., the scale will automatically indicate the new weight, you need to remain standing on the scale.
• Once you step off the scale after weighting, it will switch itself off automatically after a few seconds.

DISPLAY CODES
• Lo: Replace the battery. Type of battery CR2032 Lithium.
• oL: Overload indicator. If weight exceeds the maximum capacity of 180kg. "oL" Appears on the display screen.

CLEANING AND MAINTENANCE
• Do not immerse the appliance in water or any other liquid.
• Clean the appliance with a damp cloth. Never use harsh and abrasive cleaners, scouring pad or steel wool, which damages the device.

This product contains recyclable materials that comply with WEEE guidelines. Do not dispose of this product as unsorted waste. Please contact your local municipality for a recycling collection point.

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